

# Grain Vita Tester

---



## Application:

To measure the viability of grain ( germination capacity ). The Grain Vita Tester works fully automatically and is microprocessor controlled.

With prepared cereal samples, the measuring instrument can be used to determine the germination capacity of a variety of different cereal species in just a few minutes with the help of a reagent fluid. The Aerne Analytic Grain Cutter is used for cutting grain grains and together with the Grain Vita Tester forms an unbeatable team.

## Technical Data:

**Casing:** PVC white

**Dimensions/mm:** Depth: 280mm Width: 210 mm Height (without cover): 225 mm

**Power supply:** 230VAC

**Scope of supply:** Grain Vita Tester, 1 measuring cylinder PP 250 ml, 1 l test tube liquid (3.5 g triphenyle tetrazolium chloride) dissolved in distilled water. 1 sample holder with handle, felt and sieve. 2 plexi glass cuvettes to fill the grain halves.

## Order Details:

**Part:**

Grain Vita Test

**Order-No.**

GVT


Cuvette

GVTK

---

## Grain Cutter

---

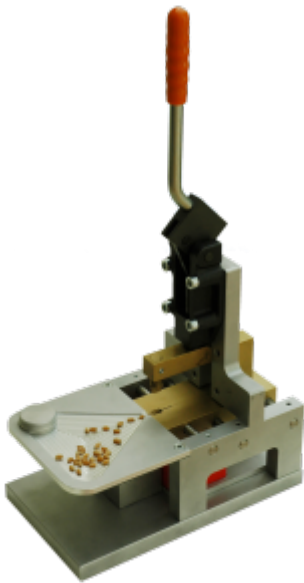

### Application:

Slicer for cutting Grain.

### Technical Data:

**Material:**

Aluminium, Brass

**Scope of supply:**

Grain Cutter with insert Plate, 2 Collecting Trays

### Order Details:

**Part:**

Grain Cutter

Bracket blade

**Order-No.**

KSN

KSNEK

---

You could also be interested in

[Mashbath](#)

